

02_01_PAA-gel-preparation (native)

FREITAG, 14.5.2021

Goal-Setting

- Preparing a PAA gel (native) to separate DNA fragments

Terms / abbreviations

- APS = Ammoniumperoxodisulfate
- PAA = Polyacrylamide
- TBE = Tris-Borate-EDTA buffer
- TEMED = Tetramethylethylenediamine

Risk areas

- Health hazard/Hazardous to the ozone layer
- Acrylamide is neurotoxic, always wear gloves!

 Hazard symbols



Required materials and / or information

- Chemicals:
 - Acrylamide 29:1; 40% (w/v), Roth
 - APS (10% (w/v))
 - MilliQ water, Sartorius arium pro VF
 - TBE Buffer (10x)
 - TEMED, AppliChem
- Materials:
 - 0.75 mm gel comb, BioRad
 - 0.75 mm gel comb, BioRad
 - 0.75 mm glass slides, one with edges and one without, BioRad
 - Falcon Tube 50 ml, Sarstedt
 - Grey sponge

Templates, devices, software

- BioRad MINI PROTEAN II Apparatus
- BioRad PROTEAN Tetra Cell

Preliminary work

- [00_02_10x-TBE-buffer-recipe](#)
- [00_01_10%-APS-recipe](#)

Operation

1. Set up the gel chamber by putting the two 0.75 mm glass slides together in the little plastic holder so that a gap between the slides is created

- Turn down the hand gears so that the glass slides are tight
- Place the two grey sponges in the big plastic holder
- Place the little plastic holders along with the glass chambers into the big plastic holder and make sure everything is really tight
 - Test it with a bit of deionized water if the chambers are leak proof
 - Dry the chambers after that with paper!
- Get all chemicals
- Dependig on the percentage of the gel, mix the components according to following:


Composition of gel for different acrylamide concentrations



	A	B	C	D
1	Component	5% PAA	15% PAA	20% PAA
2	ignore this row -->	0.05	0.15	0.2
3	40% Acrylamid [mL]	12.5	37.5	50
4	10x TBE buffer [mL]	10.00	10.00	10.00
5	H2O (TBE) [mL]	76.96	51.96	39.46
6	Mix by inverting Falcon Tube			
7	10% (w/v) APS [mL]	0.5	0.5	0.5
8	TEMED [mL]	0.04	0.04	0.04
9	total [mL]	100	100	100

- Mix the gel by inverting the Falcon Tube
- APS and TEMED will start the polymerisation so make sure that the gel chambers are ready, leakproof and dry
 - In order to get rid of air bubbles, the schottflask can be put into an UltraSound Device for a few seconds for degassing, but keep in mind that the schottflask is contaminated once it was used in the gel electrophoresis room
- Fill up the chambers to the rim of the small glass slide
- Stick the comb into the chamber and avoid bubble formation

Chamber with (left) and without (right) air bubbles

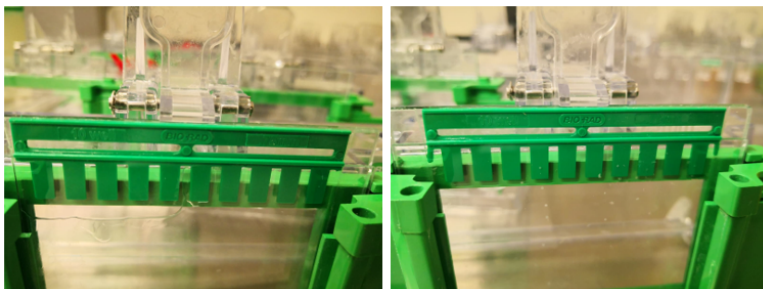

- Wait until the solution has polymerised
 - Proof the solidity with the rest of the solution in the Falcon Tube
 - This usually takse some time (~45 min)
- Take the little plastic holder out of the big one and take out the glass slides
- Gently remove the comb (not recommended if the gel should be stored)
- If the gel is not used directly, wrap it in wet paper towels als put it in a plastic bag
- Write down the date and the percentage of the gel and store at 4 °C in the fridge in electrophoresis room

Disposal

- TBE Buffer can be discarded in sink with a lot of water

Troubleshooting

- Use fresh APS
 - APS catalyzes the polymerization of acrylamide
 - Using old APS or APS stored above -20 °C will result in slow or incomplete polymerization
 - Keep small, fresh aliquots in the freezer

Follow-up work

- [02_02_Performing-PAA-gel-electrophoresis](#)