

# 03\_02\_Heat-inactivation-of-TdT-reaction

MONTAG, 27.9.2021

## Goal-Setting


- Inactivate TdT reaction through heating

## Terms / abbreviations

- TdT = Terminal deoxynucleotidyl transferase

## Risk areas

- If spilled, always wipe surface with alcohol

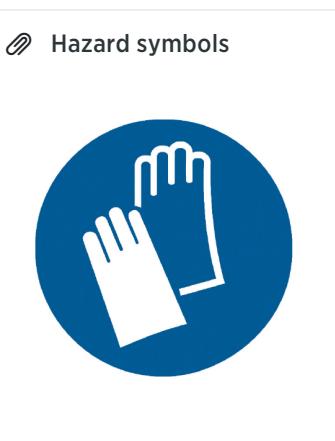

## Required materials and / or information

- Samples

## Templates, devices, software

- Thermocycler, Eppendorf


## Preliminary work

- [03\\_01\\_NEB-protocol-for-TdT-tailing-reaction](#)
- [03\\_01\\_Thermofisher-protocol-for-TdT-tailing-reaction](#)

## Operation

1. Turn on thermocycler through the switch on the back left side of the device
2. Use Benutzername and Passwort to login
3. Place samples into cycler by opening the lid. Fully close lid afterwards.
4. In iGEM21 open 'Stop' folder and start '70stop' program (70 °C, 20 min)
  - a. Always start the thermocycler in advance
  - b. Inactivating samples for 10 min should be enough
5. When finished, logout and shut down afterwards
6. Turn the electricity off when shutting down

## 🔗 Thermocycler navigation

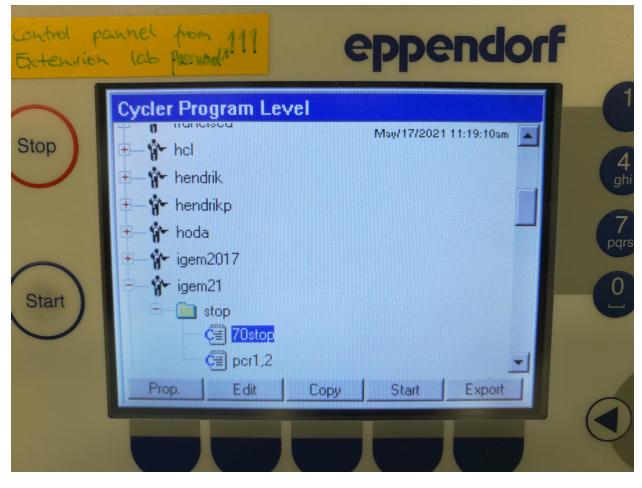

## Disposal

- None

## Troubleshooting

- If something went wrong, restart the thermocycler
- Wear gloves to reduce the risk of DNase and RNase contamination

## Follow-up work

- [01\\_02\\_Sample-preparation-for-gel-electrophoresis](#)